# NEO

## SPORT

# NS-500

## Massage gun

Pistolet do masażu

Massagepistole

**EN**

### Manual Instructions

Read and follow the instructions before installation and use. Keep the manual for later reference.

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150000 / 150001 / V. 1.2

## Warnings!

- This product is for home use only, not for commercial or medical use.
- This product is for healthy adults only, and should not be used on an injured body part. If you experience pain or discomfort, stop using the device.
- Use the massager on soft tissues of the body without causing pain or discomfort, such as muscles, tendons, ligaments and fascia.
- When the massager applies a lot of pressure or beats, bruising may occur. Check the massage area regularly and stop using the device immediately if there is inflammation, swelling or increased pain.
- Use the massager on a dry and clean area of the body by gently pressing and moving it.
- To avoid jamming, keep fingers, hair or other body parts away from the shaft or back of the massager head.
- Do not use on or around sensitive parts such as the head, face, cervical vertebrae, reproductive organs, ceramics, dentures, etc.
- Avoid excessive pressure and excessive use of this product. Otherwise, it may cause injury and shorten the life of the product.
- Inspect the device before use. If there is any damage, stop using the device. To avoid danger, do not disassemble or repair the device yourself.
- Persons with aneurysm, bleeding, cancer, heart disease or fractures, as well as infants, pregnant women and persons carrying a heart rhythm regulator or having metal contents in the body should not use this product.
- Persons with insufficient physical or mental ability, experience or knowledge, including children, should not use this product unless they do so under the supervision or direction of those responsible for their safety. Children should be supervised so that they do not play with electrical equipment.
- The surface of the device is hot and people who are not sensitive to heat should exercise caution when using it.
- Avoid direct sunlight and keep away from liquids, heat or fire.
- Do not use this product in a humid environment (such as a bathroom) to avoid the danger of short circuit.
- Choose the appropriate massage head according to your needs. Do not use the massage head on the same area for a long time.
- Do not apply excessive pressure to avoid damaging muscles and soft tissues. In addition, when excessive pressure is applied, the device will activate the safety mechanism and stop working. To continue using the device, wait 3 seconds and restart the product.
- If you experience discomfort during use, stop using the device immediately.
- Do not disassemble or modify this product in any way.
- Do not leave a running or charged product unattended.
- After 30 minutes of continuous use, turn off the product's power for 5-10 minutes to prolong its life and ensure safety.
- The product must not be used while charging.
- The power cord must not be repaired.

- Never use the device when the cord, plug or charger is damaged.
- Do not pull the power cord, do not hang the device on the power cord, avoid contact of the cord with sharp edges.
- Do not plug in or use the device with wet hands.
- Turn off the device before charging the battery.

## **Battery warnings!**

- Do not charge near a fire or in a very hot environment.
- If the charge level of the battery is low and it is not charged for a long time, the battery will go into low charge protection mode. The product will not be able to be turned on and operate normally.
- If the product is not used for a long time, it is recommended to charge it once a month to avoid damaging the battery.
- Stop using the product if the battery emits a strange odor, heats up, is deformed or shows other abnormalities.
- When the battery has reached the end of its life cycle, it should be removed from the device before disposing of it. Follow the recommendations of the waste treatment institution in your area.
- Turn off the power when removing the battery.
- Dispose of the battery in a safe manner.
- How to remove the battery: Turn and unscrew the massage head, try to open the circular decorative ring, circle the two screws, then remove the decorative stickers on the lower left and lower right corners, remove the two screws, open the case and remove the battery.

## **Charging tips**

- Before using it for the first time, charge it fully for three hours.
- When charging, connect the DC end of the USB data cable to the charging port of the battery, and then insert the cable on the AC side into the socket.
- The LED on the battery flashes regularly and displays the battery level, indicating that the battery is being charged.
- When the battery is connected, the LED on the device flashes and displays the battery level, indicating that the battery is being charged.
- When all LEDs stop flashing and are on, it indicates that the battery is fully charged.
- The battery can be charged at any time and at any level.
- The average operating time of the device is two hours, depending on the speed and the applied pressure during the use.


# Design and part list

## I Panel

- A Type c socket
- B Power indicator
- C Gear indicator
- D Switch and Intensity adjustment

## II Massage Heads

- A **U-shape:** massage of muscles near vertebrae; protection of tendons, it works on the muscles only
- B **Ball:** massage of larger muscle areas, covering a smaller area, characteristic of medium intensity
- C **Flat head:** the deepest muscle tissues with a firmly compacted structure (chest or quadriceps of the thigh)
- D **Bullet:** the most intense point massage; removing radiating muscle pains
- E **Wrist head:** massage the muscles near the bones
- F **Mushroom:** very gentle massage of the most tense parts of the muscles
- G **Banana:** massage of wide parts of the muscles
- H **Paddle:** used for larger muscle parts; smaller massage area and less intensity
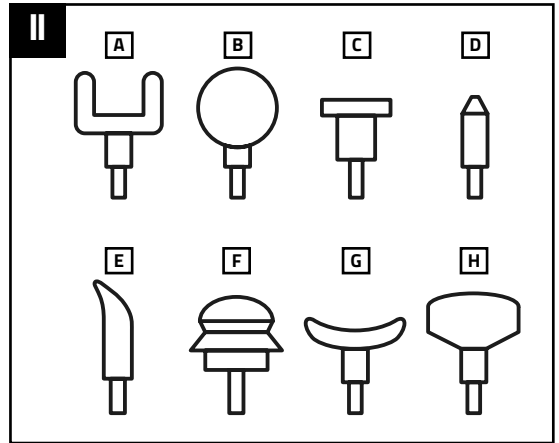
### Operation Program of Lower Limbs

#### Prone position

1. Small flat head vertically downward (The movement direction is from up to down)  
Hamstring:
  - Biceps femoris: 20s-30s × 3 times
  - Semitendinosus: 20s-30s × 3 times
  - Semimembranosus: 20s-30s × 3 times
2. Small flat head outward 45° (The movement direction is from up to down)
  - Gracilis: 10s-15s × 2 times
  - Adductor longus: 10s-15s × 2 times
  - Adductor brevis: 10s-15s × 2 times
  - Triceps surae
  - Gastrocnemius: 20s-30s × 3 times
  - Soleus: 20s-30s × 3 times
3. Small round head vertically downward  
Popliteus:
  - Knee straightening position: 15s-30s
  - Knee bending position: 15s-30s

#### Plantar fasciitis


##### Prone position (single side)

1. Small flat head vertically downward (The movement direction is from up to down.)  
Hamstring:
  - Biceps femoris: 20s-30s × 3 times
  - Semitendinosus: 20s-30s × 3 times
  - Semimembranosus: 20s-30s × 3 times
2. U head vertically downward (The front foot sole touches the ground.)  
Popliteus:
  - Knee straightening position: 15s-30s
  - Knee bending position: 15s-30sTriceps surae:

#### Supine position

1. Small flat head vertically downward (The movement direction is from up to down.)  
Quadriceps femoris:
  - Rectus femoris: 30s-45s × 2 times
  - Vastus lateralis: 30s-45s × 2 times
  - Vastus medialis: 30s-45s × 2 times
  - Vastus intermedius: 30s-45s × 2 times
2. Small flat head vertically downward
  - Sartorius: 30s-45s
  - Tensor fasciae latae: 15s-20s × 2 times
  - Tibialis anterior: 15s-20s × 2 times
  - Peroneus longus: 15s-20s × 2 times
  - Peroneus brevis: 15s-20s × 2 times

- Gastrocnemius: 20s-30s × 3 times
  - Soleus: 20s-30s × 3 times
3. Small round head vertical to planta pedis:
    - Plantar fascia
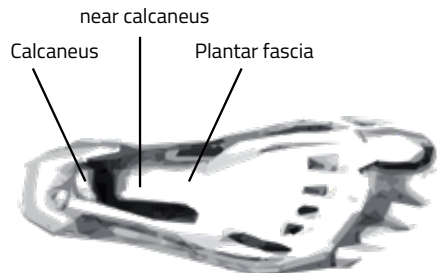
# Gonarthrititis

## Supine position

1. Small flat head vertically downward (The movement direction is from up to down and The patient slowly and slightly bends the knee.)

Quadriceps femoris:

- Rectus femoris: 30s-45s × 2 times
- Vastus lateralis: 30s-45s × 2 times
- Vastus medialis: 30s-45s × 2 times
- Vastus intermedius: 30s-45s × 2 times

2. Small round head vertically downward

- Iliac bone edge -iliac bone anterior: 30s-1min
- Rectus femoris, inguen: 1min
- Tensor fasciae latae: 15s-20s × 2 times
- Patellar ligament area: 20s × 3 times

## Prone position

1. Small round head vertically downward

- Popliteus: 30s-45s
- Knee straightening position 15s-30s
- Knee bending position: 15s-30s
- Triceps surae
- Gastrocnemius: 20s-30s × 3 times
- Soleus: 20s-30s × 3 times

# Operation Program of Sacroiliac Part

## Prone position (single side)

1. Small flat head vertically downward (The direction is from inside up to outside down, along the muscular direction.)

- Gluteus medius: 30s-1min × 4 times
- Gluteus maximus: 30s-1min × 4 times
- Piriformis: 1min-1.5min
- Articulatio sacroiliaca 1min-1.5min

2. The direction is from inside up to outside down, along the muscular direction. (From up to down)

Hamstring:

- Biceps femoris: 20s-30s × 3 times
- Semitendinosus: 20s-30s × 3 times
- Semimembranosus: 20s-30s × 3 times

## Supine position

1. Small flat head vertically downward (The movement direction is from up to down.)

- Tibialis anterior: 15s-20s × 2 times
- Peroneus longus: 15s-20s × 2 times
- Peroneus brevis: 15s-20s × 2 times

## Indications and Contraindications

The deep penetration of the device can reduce pain from joint sprains, bursitis, muscle cramps, neuralgia and other musculoskeletal disorders in a short period of time. Almost all soft tissue diseases are associated with trauma, excessive muscle use or abnormal posture, and most of the scar tissues involved in these diseases can be treated with this device. Avoid excessive stress and prolonged physiotherapy when using the device. Also avoid sensitive areas such as head, face and cervical vertebrae. Stop physiotherapy if there is inflammation, swelling or increased pain. The ongoing research is about the impact of vibration training on specific medical conditions. This is very likely to reduce the types of contraindications mentioned above, and practical experience has shown that in many cases, it is also a physiotherapy recommendation to combine vibration training into physiotherapy programs. This must be done by, or in the company of, a doctor, expert or physical therapist.

### Indications

- Pain and cramps due to muscle injury, sprains, strains
- Help the flow of edema fluid in the swollen defibrillator area
- Relax thickened connective tissues and fascias
- Reduce the accumulation of lactic acid in internal fixation.
- Increase joint mobility
- Eliminate muscle fatigue

### Contraindications

- Aneurysm, bleeding, use of blood thinner
- Heart disease, with a pacemaker or defibrillator
- Pregnancy, cancer
- Within 90 days of intra-articular fixation, the device should not be used within 3 inches of internal fixation.
- Sensitive area: head, face, cervical vertebrae, spine
- Positions close to the bones with less muscle
- tibia, the back of the foot, the back of the hand

#### **(EN) INFORMATION FOR USERS ON HOW TO HANDLE WITH WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT.**

If the appliance, packaging, operating instructions and the like are marked with the symbol of a crossed-out wheeled bin, this means that the appliance is subject to separate waste collection in accordance with Directive 2012/19/EU of the European Parliament and of the Council. The used appliance should not be disposed of with other household waste and should not be treated as municipal waste. By throwing your electro-waste into a container you pose a risk to the environment.

Take your used appliance to a separate collection point for electrical and electronic equipment organised by the public administration. By separating and submitting used electrical and electronic equipment for treatment, recovery, recycling and disposal, you protect the environment from pollution and contamination, contribute to reducing the use of natural resources and the cost of manufacturing new equipment. Proper disposal and scrapping helps to eliminate the negative impact of discarded products on the environment and human health.

For details on how this product can be recycled, please contact your local municipality, city cleaning service or the shop where you purchased the product.

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